

# ECL X2-24 ECLIPS DETACHABLE WHEELCHAIR



*For illustration purposes only*

**VERMEIREN**

## ECL X2-24 ECLIPS DETACHABLE WHEELCHAIR

### PRODUCT OVERVIEW




### ECL X2-24 Eclips Detachable Wheelchair

- Multi-function wheelchair with detachable armrests and detachable footrests for easier patient transfer.
- Breathable woven nylon seat
- Lightweight aluminium frame

*For illustration purposes only*

# DETACHABLE WHEELCHAIR

## PRODUCT FEATURES

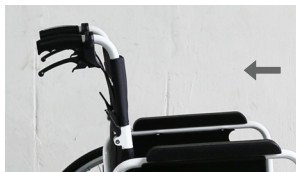

### FOLDABLE BACKREST

Press on the lever to release and fold back rest downwards for both sides



### DETACHABLE ARM REST

Press the release button to unlock and remove the armrest.



### DETACHABLE FOOTREST

Pull the lever to release the lock and remove the footrest.



### STATIONARY HAND BRAKE

Push up to unlock and down to lock

*For illustration purposes only*

## CARE INSTRUCTIONS

- When using the Wheelchair, please pay attention to people and/or items around to avoid collision.
- Lock the wheels before getting in and out of the Wheelchair.
- When getting in and out of the Wheelchair, please do not stand on the footplates.
- The footplates should be folded up / detached beforehand.
- Pay attention to the shift of the Center of Gravity when on elevated grounds.
- When passing through doors or arches, make sure that there is enough room to avoid causing injuries to the user, or the Wheelchair.
- Do not exceed the maximum safe working load.

## ECL X2-24 ECLIPS DETACHABLE WHEELCHAIR

### PRODUCT SPECIFICATIONS

|                                     |                |
|-------------------------------------|----------------|
| Overall Length                      | 99 cm or 39"   |
| Overall Width (Wheel to Wheel)      | 66 cm or 26"   |
| Overall Height                      | 92 cm or 36.2" |
| Folded Length                       | 80 cm or 31.5" |
| Folded Width                        | 37 cm or 14.6" |
| Folded Height                       | 76 cm or 29.9" |
| Seat Depth                          | 41 cm or 16.1" |
| Seat Width                          | 45 cm or 17.7" |
| Seat Height                         | 46 cm or 18.1" |
| Arm Rest Height (Floor to Arm Rest) | 67 cm or 26.4" |
| Back Rest Height                    | 46 cm or 18.1" |
| Front Wheel Diameter                | 15 cm or 5.9"  |
| Back Wheel Diameter                 | 59 cm or 23.2" |
| Safe Working Load                   | 120 kg         |
| Overall Weight                      | 15 kg          |

Check out our full range of detachable wheelchairs. When comparing models, consider the user's environment and how tasks such as transfers will be carried out before deciding

Scan the QR code to find out more about the product. To scan, use your camera app or a QR code reader on your device.

