

FT6802 AIRPLANE WHEELCHAIR



FT6802 AIRPLANE WHEELCHAIR

PRODUCT OVERVIEW



Lightweight Wheelchair

- Easy to fold and operate
- Multi-function wheelchair with fold-up foot rests for easier patient transfer.
- Foldable frame for easy storage
- Ventilated and cushioned seat upholstery
- Aluminium lightweight frame
- Pull-Trolley attachment for easy mobility


LIGHTWEIGHT WHEELCHAIR

PRODUCT FEATURES



FLIP UP ARM REST

Flip up arm rest by pushing up

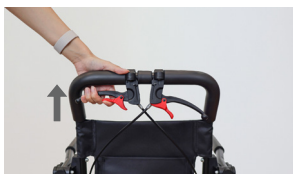

HANDLE BRAKES

When brakes are engaged, press the red trigger upward to release and move the wheelchair.



FLIP UP FOOT REST

Hold the seat frame and backrest tube. Gently pull inward to close.

For illustration purposes only

CARE INSTRUCTIONS

- When using the wheelchair, please pay attention to people and/or items around to avoid collision.
- Lock the wheels before getting in and out of the Wheelchair.
- When getting in and out of the Wheelchair, please do not stand on the footplates.
- The footplates should be folded up / detached beforehand.
- Pay attention to the shift of the Center of Gravity when on elevated grounds.
- When passing through doors or arches, make sure that there is enough room to avoid causing injuries to the user, or the Wheelchair.
- Do not exceed the maximum safe working load.

FT6802 AIRPLANE WHEELCHAIR

PRODUCT SPECIFICATIONS

Overall Length	88 cm / 34.6"
Overall Width (Wheel to Wheel)	51.5 cm / 20.3"
Overall Height	93 cm / 36.6"
Folded Length	54.5 cm / 21.5"
Folded Width (Wheel to Wheel)	43 cm / 16.9"
Folded Height	71 cm / 28"
Seat Depth	36 cm / 14.2"
Seat Width	40.5 cm / 15.9"
Seat Height (Floor to Seat Base)	51 cm / 20.1"
Arm Rest Height (Floor to Arm Rest)	73 cm / 28.7"
Back Rest Height	34 cm / 13.4"
Front Wheel Diameter	15 cm / 5.9"
Back Wheel Diameter	30.5 cm / 12"
Safe Working Load	80 kg
Overall Weight	7.6 kg

Check out our full range of travel wheelchairs. When comparing models, consider the user's environment and how tasks such as transfers will be carried out before deciding

Scan the QR code to find out more about the product. To scan, use your camera app or a QR code reader on your device.

