

# PRACTICO 25/80 PREMIUM FOLDABLE RAILS NURSING BED




# PRACTICO 25/80 PREMIUM NURSING BED, FOLDABLE SIDE RAILS

## PRODUCT OVERVIEW

Hospital Bed Purchase Includes:

- 4" Premium Comfort Mattress

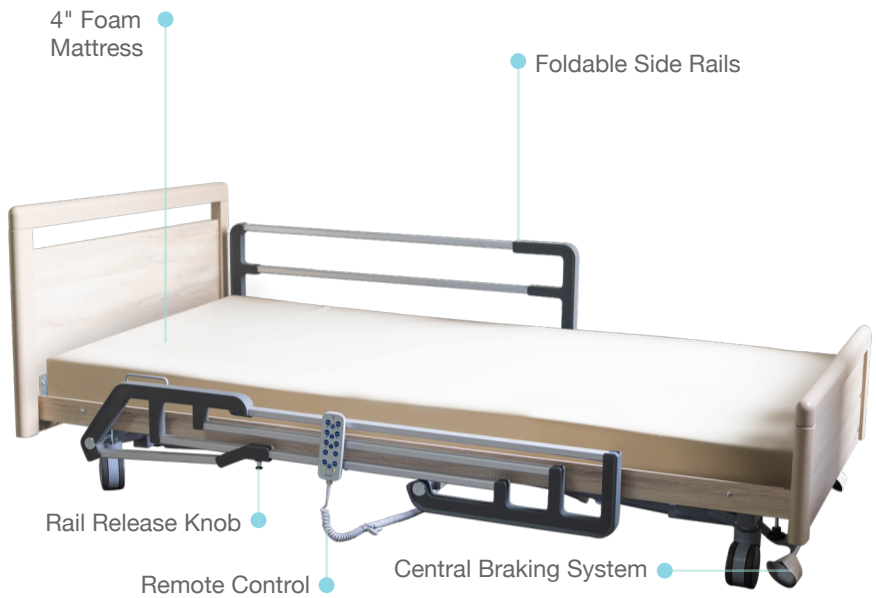

*For illustration purposes only*



Foldable Side Rail



Rail Release Knob



Central Braking System




Remote Control

*For illustration purposes only*

## PRACTICO 25/80 PREMIUM NURSING BED, FOLDABLE SIDE RAILS

### PRODUCT FEATURES

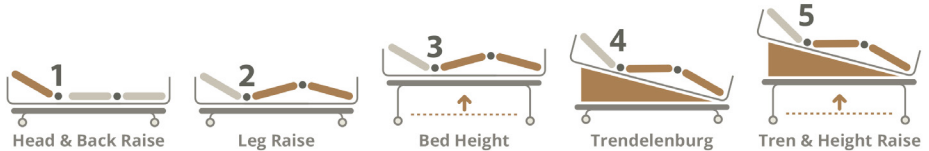

The 5 Crank electrical premium bed is our standard electrical bed for home patient use. It is low in power consumption and only consumes power when the remote is pressed on.

#### Bed Functions:

- Back Raise lifting angle is 0 - 70° ( $\pm 5^\circ$ )
- Foot Raise lifting angle is 0 - 50° ( $\pm 5^\circ$ )
- Trendelenburg Position is 0 - 15° ( $\pm 5^\circ$ )



All steel parts are powder coated. Each bed can be equipped on both side with side rails. Side rails can be lifted and lowered to protect users against falling out of bed, at the same time enable safe and comfortable boarding and alighting.

## CARE INSTRUCTIONS

- Ensure adequate space around the bed to prevent user injury or bed damage
- Prior to entering or exiting the bed, confirm that the wheels are securely locked
- Seek assistance when adjusting the bed if users encounter difficulties
- Do not surpass the specified maximum weight capacity

## PRODUCT SPECIFICATIONS

### PRACTICO 25/80 BED, FOLDABLE SIDE RAILS

Key Features	Head/Back Raise • Leg Raise • Adjustable Height Trendelenburg • Trendelenburg & Height Adjustment
Bed Railings	Foldable Side Rails
Overall Length	210 cm / 82.6"
Overall Width	105 cm / 41.3"
Overall Height	25 - 80 cm / 9.8" - 31.4" (Floor to Bed base)
Head/Back Raise Angle	0° to 70°
Leg Raise Angle	0° to 50°
Trendelenburg Position	0° to 15°
Safe Working Load	185 kg
Country of Origin	Made in Germany

Check out our full range of Hospital Beds & Mattresses. When making comparisons, our best suggestion would be for you to think about your environment and consider how certain task like transferring someone is done before deciding.

Scan the QR code to find out more about the product. To scan, use your camera app or a QR code reader on your device.

