WALKING AIDS | USER GUIDE

rainbow care

FS965LH-B ROLLATOR



For illustration purposes only

rainbow care

FS965LH-B ROLLATOR

PRODUCT OVERVIEW



PRODUCT OVERVIEW



For illustration purposes only

FS965LH-B Rollator

The FS965LH-B rollator has a flip-up seat and footrest, allowing the user to rest when needed. Ensure the brakes are locked before sitting down. The basket provides convenient storage.

rainbow care

FS965LH-B ROLLATOR

PRODUCT FEATURES






REMOVE THE BASKET


Lift up the seat and pull up the basket



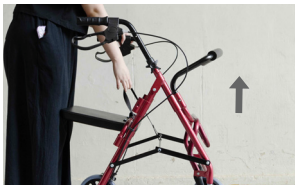

FOLDABLE

Remove the basket and pull the string in the middle upward to fold

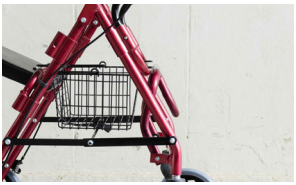



FLIP UP FOOT REST

Flip up to keep or flip down for usage





CHANGEABLE SEAT SUPPORT

Press the button and pull it up to change the direction

For illustration purposes only

CARE INSTRUCTIONS

- When using rollator, please pay attention to people and/or items around to avoid collision.
- Make sure the level of the height of each handle is the same before using it.
- Make sure the level of flooring is even before seating down.
- Do not exceed the maximum load.

rainbow care

FS965LH-B ROLLATOR

PRODUCT SPECIFICATIONS

Overall Length	78.5 cm / 30.9”
Overall Width	58.5 cm / 23”
Overall Height	77 - 88 cm / 30.3 - 34.6”
Folded Length	34.5 cm / 13.6”
Folded Width	58.5 cm / 23”
Folded Height	90.5 cm / 35.6”
Seat Depth	31.5 cm / 12.4”
Seat Width	34.5 cm / 13.6” (Cushion)
Seat Height	53.5 cm / 21.1” (Floor to Seat)
Front Wheel Diameter	20 cm / 7.9”
Back Wheel Diameter	20 cm / 7.9”
Safe Working Load	80 kg 5 kg (Basket)
Overall Weight	8.8 kg

Check out our full range of rollators. When comparing models, consider the user’s environment, mobility needs and required safe working load.

Scan the QR code to find out more about the product. To scan, use your camera app or a QR code reader on your device.

6223 0904 • 6931 7425 8501 1878 sales@rainbowcare.com.sg rainbowcare.com.sg
@rainbowcare rainbowcaresg Blk 3007 Ubi Road 1 #02-410 Singapore 408701