WALKING AIDS | USER GUIDE

rainbow care

FS966LH-B ROLLATOR



For illustration purposes only

rainbow care

FS966LH-B ROLLATOR

PRODUCT OVERVIEW



PRODUCT OVERVIEW



For illustration purposes only

FS966LH-B Rollator

The FS966LH-B rollator has a flip-up seat and footrest, allowing the user to rest when needed. Ensure the brakes are locked before sitting down.

rainbow care

FS966LH-B ROLLATOR

PRODUCT FEATURES





REMOVE THE BASKET


Lift up the seat and pull up the basket





FOLDABLE


Remove the basket and pull the string in the middle upward to fold

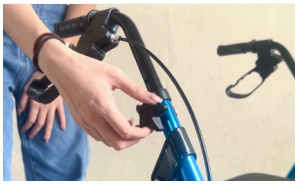



HEIGHT ADJUSTMENT

Turn the knob, remove the screw to change the backrest

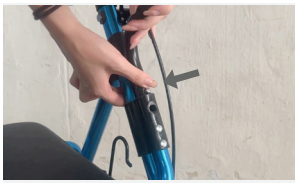



CHANGEABLE SEAT SUPPORT

Press the button and pull it up to change the direction

For illustration purposes only

CARE INSTRUCTIONS

- When using rollator, please pay attention to people and/or items around to avoid collision.
- Make sure the level of the height of each handle is the same before using it.
- Make sure the level of flooring is even before seating down.
- Do not exceed the maximum load.

rainbow care

FS966LH-B ROLLATOR

PRODUCT SPECIFICATIONS

Overall Length	63 cm / 24.8”
Overall Width	60 cm / 23.6”
Overall Height	79.5 - 89 cm / 31.3 - 35”
Folded Length	49 cm / 19.3”
Folded Width	60 cm / 23.6”
Folded Height	91 cm / 35.8”
Seat Depth	32.3 cm / 12.7”
Seat Width	35.5 cm / 14” (Cushion)
Seat Height	50.5 cm / 19.9” (Floor to Seat)
Front Wheel Diameter	15 cm / 6”
Back Wheel Diameter	15 cm / 6”
Safe Working Load	80 kg 5 kg (Basket)
Overall Weight	8.5 kg

Check out our full range of rollators. When comparing models, consider the user’s environment, mobility needs and required safe working load.

Scan the QR code to find out more about the product. To scan, use your camera app or a QR code reader on your device.

6223 0904 • 6931 7425 8501 1878 sales@rainbowcare.com.sg rainbowcare.com.sg
@rainbowcare rainbowcaresg Blk 3007 Ubi Road 1 #02-410 Singapore 408701